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The after services of three guarantees (for repair, replacement and compensation of faulty products) of DB-Phototherapy apparatus is strictly in accordance with the Law of the People's Republic of China on the Protection of Consumer Rights and Interests and the Law of the People's Republic of China on Product Quality. The service contents are as follows:

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### Warranty ordinance:

- 1.The warranty period is 12 months from the date of purchase.
  - 2.Within 7 days from the day after you sign for the product, if the product has any performance failure, you can enjoy free return or replacement services after the product is tested and confirmed by DB-light Technology after-sales service center.
  - 3.Within 8-15 days from the day after you sign for the product, if the product has any performance failure, you can enjoy free replacement or maintenance services after the product is tested and confirmed by DB-light Technology after-sales service center.
  - 4.Within 12 months from the day after you sign for the product, if the product has any performance failure, you can enjoy free maintenance services after the product is tested and confirmed by DB-light Technology after-sales service center.
  5. You can enjoy free maintenance services if there was any material defect, malfunction or technological problem in the using process following to the instructions during the warranty period.
  6. There will be a maintenance cost if the following situation occur within the warranty period:
    - Product damage caused by flood, fire disaster, abnormal voltage and other natural disaster.
    - Faulty or damage caused by using wrongly and Unauthorized maintenance, renovation.
    - Can not present the warranty card and proof of purchase 365 days after the date of production.
    - The user information of the warranty card is not clearly or correctly recorded.
    - Product failure caused by falling off or transportation.
    - Product damage caused by Unauthorized maintenance or disassembly without the testing by the company's personnel test sales service personnel or the company's technical maintenance personnel.
- Non-warranty ordinance:
- 1.Unauthorized maintenance, modification, collision; negligent, improper use; misuse non-product accessories.
  - 2.The period of three guarantees has expired.
  - 3.Damage caused by force majeure.
  - 4.Performance failure of this product and its accessories caused by human reason.

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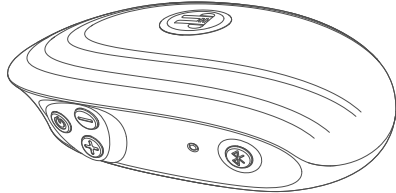
## ● Company profile

Zhongshan DB-light Technology Co.,Ltd. is the first high-tech enterprise incubated by Zhongshan-Fudan Joint Innovation Center founded by Zhongshan City and Fudan University jointly, that mainly engage in optical health and optical medical. DB-light Technology Co.,Ltd. integrates production, scientific research, production, sales and operation service system, serves the purpose of human health.

Zhongshan DB-light Technology Co.,Ltd. relies on the strong research & development capabilities of Fudan University and the industrialization of a large number of scientific and technological inventions and patented technologies to service in the cause of human health.

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## ● Technical introduction

Innovatively breaking through the theory of 'Biphasic Effect' by professor Hamblin of Harvard University and innovating cellular optical deep repair technique, precisely controlling the light and taking it into myocytes deeply.

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## ● Technical features

Deeply penetrating directly to the myocytes and realizing the regulation and repair of cell growth by light-sensitive effects of CCO and other substances in mitochondria and the effects on ATP, MMP, ROS and NO in cells. It is an intelligent and movable belly button phototherapy instrument with biophotochemical change process. Its unique design makes it portable and easy to use.

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## ● Product value

Effectively relieve navel cold, abdominal discomfort, endogenous dampness, endogenous cold, etc.

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## ● Preliminary instructions

Thank you for buying our product! Please be sure to use the product correctly according to the instructions and suggestions.

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### Caution:

- 1.This product is a health care device but not a medical device.
- 2.Do not look at the ray of light directly for a long time in case of affecting your eyesight.
- 3.Check whether the skin at the phototherapy site has normal sense of temperature, and strict observation is recommended for user with disorders.
- 4.If there is wound on the irradiated site, clean first.
- 5.You may sweat a lot during use. It is recommended to drink 300-500ml of water 30-60 minutes before phototherapy.
- 6.Red light irradiation should not be used locally within 24-48 hours where an acute trauma occurs, so as not to aggravate swelling, pain and bleeding. Low-intensity exposure is feasible after the acute phase.
- 7.Do not use the product if it is damaged. Contact the manufacturer WeChat customer service or public customer service in time, and upload pictures, communicate in time to solve the problem.
- 8.Radio waves generated by remote control may affect the normal operation of implantable medical devices or personal medical devices, such as pacemakers, cochlear implants, hearing AIDS, etc. If you use any of these medical devices, please consult their manufacturer about the limitations of using these devices.
- 9.It is recommended to use the product within the temperature range of -10℃~+40℃ and to store the product and its accessories within the temperature range of -40℃~+70℃. Equipment failure may occur if the ambient temperature is too high or too low.

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## Safety Caution:

- Do not use this product if there was any damage on the lens or the shell.
  - Do not dismantle, repair or remould it unauthorizodly.
  - Do not use it on children and keep it out of their reach.
  - Please remove the USB charging cable from the socket when it is not used.
  - Stop using it and shut it down if there was something wrong (like too hot, the smell of burning).
  - If you have an anaphylactic reaction when using it, stop using immediately and consult a doctor.
- The following persons are not allowed to use this product:
- 1.Intemperant and person who is drunk.
  - 2.The patient with xeroderma pigmentosum(DNA repair enzyme deficiency).
  - 3.The patient with unstable blood pressure(or blood pressure goes up and down frequently) or abnormalities in heart rate.
  - 4.The patient with acute trauma within 24 or 48 hours.
  - 5.Not interded for use by pregnant or nursing women.

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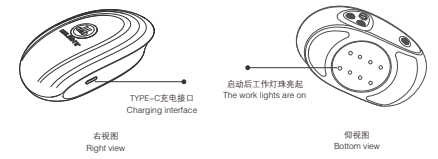
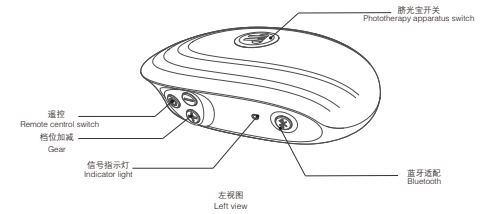
蓝牙连接用于控制本公司其它设备  
Bluetooth remote control function is used to control other equipment of the company.

 脐光宝光疗仪可以针对身体肚脐部位使用  
对肠胃等均有显著效果  
Phototherapy apparatus can be used for navel of the body.  
Phototherapy have a significant effect on gastrointestinal care.

透光区贴近皮肤启动光疗  
Start phototherapy by pressing the light zone close to the skin.

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#### Common problems in the course of using

##### 1.How to use the product?

Clean and dry your skin is recommended before using the product. In order to guarantee the use effect,there should be no clothes between the product and your skin. The soft surface of the product is made of high quality material, which can contact with the skin directly.

##### 2.Is there any requirements about the frequency of using the product?

Using once per day for 10 minutes at a time is recommended. The product will automatically shut down when the single use time is reached.

##### 3.The cautions during the use of the product.

- If there is a warm feeling on your skin during use, it is normal while the produce is working.
- Do not look at the ray of light directly for a long time in case of affecting your eyesight.

##### 4.What if the product and essential oil are used together?

Colorless essential oil is recommended if you want. The essential oil with color will affect the phototherapy effect of the product like light clothes.

##### 5.What if the product is cleaned with water?

The product cannot be washed by water. It's shell is made of high quality material, you can wipe it with soft cloth if there are any blot on it. If the product fall into the water by accident, take it out immediately in case of a short circuit causing fire.

##### 6.what is the efficacy of the red light?

According to the research about Low level light therapy(LLLT) and photobiomodulation(PBM), the irradiation of red light on deep tissue can reduce inflammation that causes pain in the cervical spine and suppress and eliminate pain triggers by make a beneficial effects on Cellular oxygen consumption, ATP produce, cell membrane potential, ROS and many other key factors.

Aiming at the problem of biphasic effect(there is no effect until the light intensity reached a certain level. The effect is poor if the light intensity is weak, but the effect will also be affected if the light intensity is too strong.) existing in traditional phototherapy equipment, a unique regulative way of light is used in DB-Phototherapy apparatus, so that the light intensity can be controlled precisely and there is no biphasic effect any more. The phototherapy effect was significantly improved in this way.

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#### Maintenance and care of the product after using:

- Long press on the power button ‘DB’ to turn it off.
- Wipe the lens with a soft cloth, as much as possible to avoid excessive force cause damage to the lens.
- Please place it properly after use.

#### Storage instructions:

- The product is not waterproof, keep it away from the bathtub or sink.
- Place the product where children can not reach.

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#### Fault treat suggestion:

Fault phenomenon	Possible cause	Solution
Can not turn on	The battery is low	Recharge
The side button remote control operation no respond	This function does not control this device	Use after adapt with other Bluetooth remote control devices of our company

#### If there are any other problems not on the list, try the following steps:

- Long press on the power button ‘DB’ until the product shut down.
- Keep the product off for 60 minutes.
- Restart the product.

If the problems are still there, please contact our customer service center for help.

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- Phototherapy apparatus
- TYPE-C charging cable
- Specification

#### Product specifications:

1.technical specifications  
Product name: QGB Phototherapy apparatus  
Product model: DBLIGHT-QGB01  
Rated voltage: 5V  
Rated power: 3W  
Rated Current: 0.6A  
Time rating: 10mins  
Battery capacity: 800mAh  
Power supply mode: polymer lithium battery  
Charging time: about 80 mins  
Charging Power: 10W(SV/2A)  
Net weight: 80g  
Product size: L:104mm、W:66mm、H:32mm  
Package size : L:151mm、W:106mm、H:42mm  
Optical wavelength: 650-670nm

#### 2.Transportation and storage conditions

Ambient temperature range: -40℃~+70℃  
Relative humidity range: ≤95%

Keep in a cry, ventilated place.

#### Executive standard:

GB/T 26572  
GB4706.10-2008  
IEC/EN 62471、CIE S009：2002、GB/T 20145  
GWJ015-2016

#### Listing of environmental protection:

- Do not treat the product and its accessories as ordinary household garbage.
- Please follow the local laws for the disposal of this product and its accessories, and support electronic product garbage collection classification.

#### The name and content of hazardous substances in the product

Parts	Hazardous substances					
	Plumbum(Pb)	Cadmium(Cd)	Mercury(Hg)	Hexavalent chromium(Cr(VI))	Polybrominated biphenyls(PBB)	Poly brominated diphenyl ethers(PBDE)
PCBA	X	O	O	O	O	O
Charging line	O	O	O	O	O	O
Shell&Components	O	O	O	O	O	O

#### This form is based on SJ/T 11364 specification

O: means that the content of the hazardous substance in all homogeneous materials of the component is below the limit specified in GB/T 26572.

X: means that the content of the hazardous substance in at least one homogeneous material of the component exceeds the limit specified in GB/T 26572, and there is no mature alternative in the industry at present, while it meets the environmental protection requirements of EU RoHS directive.

#### Disclaimer:

Zhongshan DB-light Technology Co.,Ltd. reserves all the right for the product explanation. If any third party legal person uses the information related to this product without the authorization of the Company, the Company will file a legal lawsuit and claim for loss caused to the Company based on such behavior. Zhongshan DB-light Technology Co.,Ltd. shall not assume any responsibility for any wrong use by any user under any circumstances.

#### Company information:

Zhongshan DB-light Technology Co.,Ltd.  
Address: 14th Floor, South Wing, Shumao Building, Torch Development Zone, Zhongshan City, Guangdong Province  
Service tel: 400 088 5099  
Landline: (86) 0760-8992 0870  
Official website: www.dblight-tech.com